A E N O

Quick Start Guide


Dental Irrigator D11




(Model: ADI0001)

v 3.1.4 | 10.09.2026


ENG.....	8
ARA.....	13
BOS.....	17
BUL.....	22
CES.....	28
DEU.....	33
ELL.....	39
EST.....	44
FRA.....	49
HRV.....	55
HUN.....	60
HYE.....	65
ITA.....	70
KAT.....	75
KAZ.....	80
LAV.....	85
LIT.....	90
NLD.....	95
POL.....	101
POR.....	107
RON.....	113
RUS.....	119
SLK.....	125
SLV.....	131
SPA.....	136
SRP.....	142
UKR.....	148
UZB.....	154

A

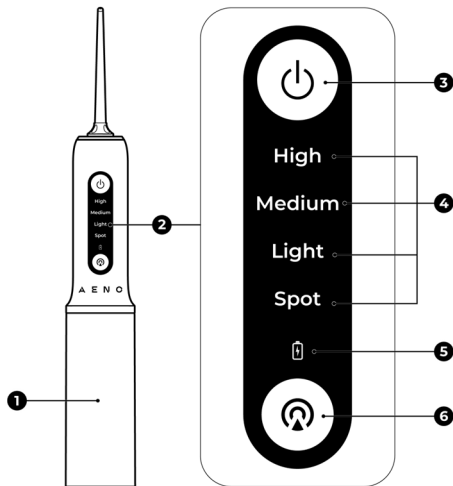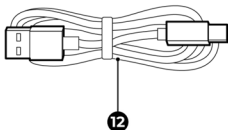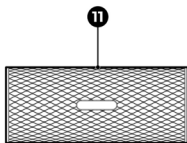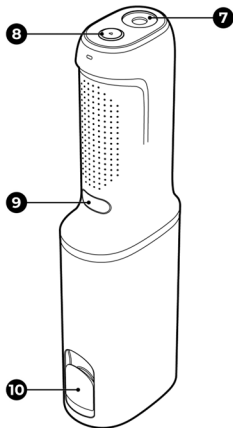C

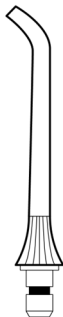
13


14


15

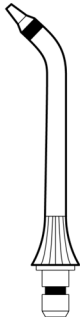
16

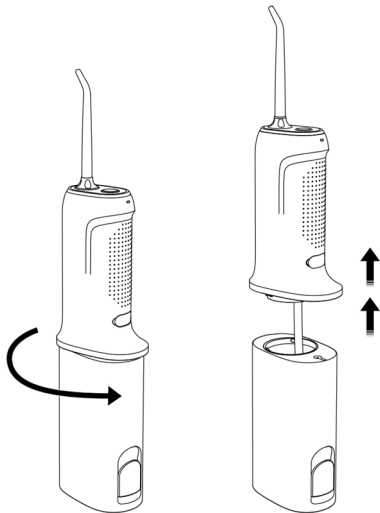E

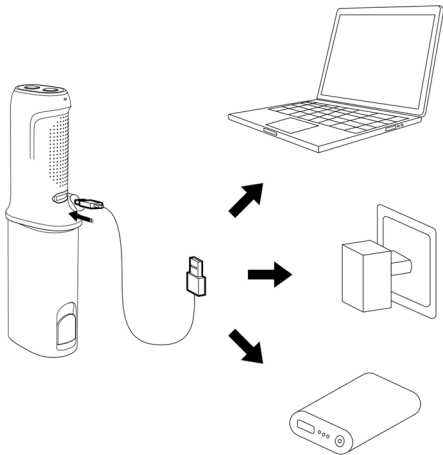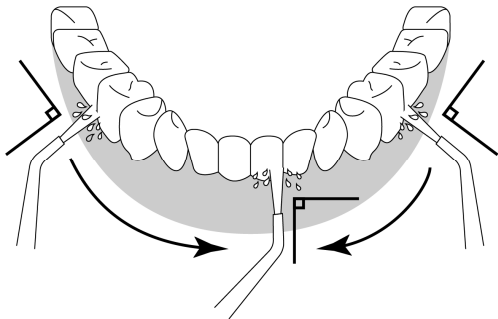ENG ADI0001 Dental Irrigator is designed for cleaning plaque and food debris in hard-to-reach places of the oral cavity.

Technical Specifications

Power supply (input): 5.0 V / 1.0 A (DC). Battery: lithium-ion, non-removable, 3.7 V, 2,200 mAh. Water pressure: 206–827 kPa. Pulsation: 1,500 pulsations/min. Full battery charge time: up to 5 h. Water tank: 200 mL. Ingress protection rating: IPX7. Operating time on a single charge: no less than 45 min. Operating time with a full tank: up to 2 min. Control: manual. USB (A) / USB (C) cable length: 1 m. Color: White. Casing material: ABS plastic. Dimensions (L×W×H): 70×45×223 mm. Noise level (at a distance of 0.3 m): ≤ 69 dB. Operation: temp. +5...+40 °C, RH 45–75 %. Storage: temp. 0...+50 °C, RH 45–90 %.

Scope of Supply (see Fig. **B, C**)

ADI0001 Dental Irrigator, replacement nozzles (4 pcs.), USB (A) / USB (C) cable, case for replacement nozzles, Quick Start Guide.

Device Elements and Accessories (see Fig. **A, B, C**)

1 – water tank, 2 – control panel, 3 – power button, 4 – operating modes indicators, 5 – charge level indicator, 6 – mode and strength selection button, 7 – opening for the nozzle, 8 – nozzle release button, 9 – USB (C) port with a cover, 10 – water tank cap, 11 – case for replacement nozzles, 12 – USB (A) / USB (C) cable, 13 – standard nozzle, 14 – tongue cleaner nozzle, 15 – orthodontic nozzle, 16 – periodontal nozzle.

Limitations and Warnings

The device is not intended for use by children or by persons with physical, sensory, or mental disabilities unless they have sufficient experience or knowledge of how to use the device, and are supervised by a person responsible for their safety. Do not allow children to play with the device and its accessories. If you have oral cavity diseases, you should consult your dentist before using the device.

If you have any questions or difficulties in using your AENO™ device, please e-mail support at **support@aeno.com** or chat with us online at **aeno.com/service-and-warranty**. The specialists will help you solve the problem, and you will not have to waste time and effort visiting the store.

Do not use a damaged device and its accessories. If the device has been transported or stored at low temperatures, leave it at room temperature for 2 hours before use. Do not drop or throw the device. Do not clean the device until it has been turned off and disconnected from the power supply. Do not use chemical and aggressive detergents, abrasive pastes, products containing acids and solvents or metal sponges for cleaning. Do not wash the device or its accessories in a dishwasher. Use original accessories only. Replace the nozzles at least once every 3 months. Do not point the outlet of the nozzle into your eyes while the device is operating. If you are not using the device for a long time, charge it every 2–3 months to prevent deep discharge of the battery. A detailed description of the device, its modes and functions can be found in the complete Operating Manual available at aeno.com/documents.

Intended Purpose of Nozzles

- The standard nozzle (13) is designed to remove plaque and food debris from interdental spaces;
- The tongue cleaner nozzle (14) is designed to remove plaque from the tongue surface;
- The orthodontic nozzle (15) is designed to clean braces;
- The periodontal nozzle (16) is designed to clean and rinse periodontal pockets.

Functions of the Device Buttons

Name	Action	Meaning
Power button (3)	Press once	Turn the device on/off
	Press and hold	Water supply in the “Spot” mode
Mode and strength selection button (6)	Press once	Select the operating mode of the device
	Press and hold for 3 seconds	Move to the next strength level
Nozzle release button (8)	Press and hold	Unlock the nozzle to remove it

Battery Charge Indication

Color	Status	Meaning
Green	On	The device is connected to a power supply and is fully charged
Red	Blinking slowly *	The device is connected to a power supply. Charging is in progress
		The device is turned on. Charge level is at 5-10 %
	Blinking fast **	Charge level is at < 5 %, the device does not turn on
–	Off	The device is turned off or not connected to a power supply
		The device is turned on. Charge level is at > 10 %

* Blinking slowly – the indicator (5) lights up for 2 seconds, then goes out for 0.5 seconds.

** Blinking fast – the indicator lights up for 0.5 seconds, then goes out for 0.5 seconds.

Operating Modes of the Device

- “High” – intensive cleaning. Recommended for advanced users;
- “Medium” – daily standard cleaning;
- “Light” – gentle cleaning for sensitive gums. It is also recommended for those who are using the device for the first time;
- “Spot” – thorough treatment of the oral cavity, removal of food debris from interdental spaces, prevention of tartar formation. Press and hold the power button for water supply in this mode.

Each operating mode has three strength levels: maximum, normal, minimum. The mode and its strength level, set during last use, will be automatically selected the next time the irrigator is turned on.

Preparing the Device for Use

1. Carefully open the packaging and remove the device and its accessories.
2. Remove any packaging materials. Remove the protective film from the control panel (2).
3. Disconnect the tank (1) as shown in Figure D. Rinse the tank and nozzles under running water, then wipe them with a soft, dry cloth.

4. Place the device on a flat, stable surface. Open the cover of the USB (C) port (9) and connect the device to a power supply with the USB (A) / USB (C) cable (12) (see Fig. E).

Note. The power adapter is not included in the package. Output parameters of the adapter must be 5.0 V / 1.0 A (DC).

5. Fully charge the device (at least 5 hours). When the battery is fully charged, the charge level indicator will turn green.
6. Disconnect the device from the power supply and insert the cover into the USB (C) port of the device.

Use

1. Open the tank cap (10) and fill the tank with water.

ATTENTION! Use purified water only without mouthwashes and other liquids or substances.

2. Close the tank cap.
3. Select the required nozzle (see "Intended Purpose of Nozzles") and insert it into the opening (7) until it clicks into place.

Note. The nozzles can rotate 360°.

4. Rotate the outlet of the nozzle away from yourself and turn the device on by pressing the power button once.

Note. The first time you turn on the device, it will operate in "High" mode at a maximum strength level.

5. Select the required mode and its strength with the mode and strength selection button (see "Operating Modes of the Device" and "Functions of the Device Buttons"). The name of the selected mode will be displayed on the control panel of the device.

Note. There is no strength level indicator on the device. The current strength level can only be distinguished during operation by the head of the supplied jet of water.

6. After selecting the required operating mode and strength, turn off the device by pressing the power button once.

7. Point the outlet of the nozzle into the oral cavity and turn on the device by pressing the power button again.
8. Perform oral hygiene by moving the outlet of the nozzle along the gum contour as shown in Figure F.

Note. If the tank runs out of water during operation, turn off the device by pressing the power button once and refill the tank as described in items 1 and 2 above. Then continue with your oral hygiene.

9. Turn off the device by pressing the power button once.
10. Open the tank cap and drain the remaining water.
11. Clean the nozzle and tank (see "Daily Care").

Cleaning and Maintenance

Daily Care. Press and hold the nozzle release button and pull the nozzle upwards. Detach the tank as shown in Figure D. Rinse the nozzle and tank under running water, then dry them thoroughly. Reinstall the tank and nozzle.

Total Care. Perform the total cleaning of the water tank once a week. Detach the tank as shown in Figure D. Wash it outside and inside under running water with a soft brush or sponge, without detergents or cleaning agents, and dry thoroughly. Reinstall the tank.

Troubleshooting

The device does not turn on. Possible cause: the battery is low. Solution: charge the device (charging time – at least 5 hours).

The device does not supply water. Possible causes: the tank is out of water; the outlet of the nozzle is clogged. Solution: fill the tank with water; clean the outlet of the nozzle.

Insufficient water supply. Possible cause: the wrong operating mode / strength has been selected. Solution: set the required operating mode / strength with the mode and strength selection button.

ATTENTION! If none of the suggested methods have solved the problem, contact your supplier or an authorized service center. Do not disassemble the device or attempt to repair it yourself.

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ENG Information contained herein may be changed without prior notice to users. For the latest information about the device, certificates, warranty, and quality, see the complete User Manual at aeno.com/documents. This device complies with applicable safety requirements.

Warranty. The warranty period and service life of the device are 2 years from the date of retail sale of the device, unless otherwise required by local law. Additional warranty information is available at aeno.com/service-and-warranty.

Manufacturer: ASBISC Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Tel.: +48732080077, info@asbis.com, asbis.com. Made in China. The date of manufacture is indicated on the package. All trademarks and trade names used herein are the property of their respective owners.

RoHS Directive. The device complies with the requirements of the Directive 2011/65/EU on the Restriction of Hazardous Substances (RoHS), including the requirements of the RoHS Directive 2015/863/EU.

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ENG Disposal information. These symbols indicate that when disposing of the device, its batteries and accumulators, as well as its electrical and electronic accessories, you must follow the rules for handling waste electrical and electronic equipment (WEEE) and the rules for handling battery and accumulator waste. According to the regulations, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories with unsorted household waste, as this will harm the environment. This device has a non-removable battery; at the end of the battery's service life, the device must be disposed of. To dispose of this device, return it to the point of purchase or take it to a local recycling center. For detailed information, contact your local household waste disposal service.

For devices with non-removable batteries: under normal or abnormal operating conditions, as well as in the event of a single malfunction, the charging voltage and current must not exceed their maximum permissible values. If the battery temperature exceeds the maximum allowed temperature, the device will stop charging under abnormal conditions. If the battery temperature drops below the minimum allowed temperature during charging, the current will be limited under abnormal conditions.

In case of battery deformation or leakage, immediately stop using the device. Put on latex/rubber gloves and a protective mask and open the windows in the room. Place the device with the leaking battery in one or two transparent plastic bags. Dispose of the device in accordance with local recycling regulations.

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