A E N O

Quick Start Guide


Sonic **Toothbrush** DB9M/DB9V


(Model: ADB0009)

v 1.0.1 | 02.02.2026


ENG.....	7
ARA.....	12
BOS.....	16
BUL.....	21
CES.....	26
DEU.....	31
ELL.....	36
EST.....	41
FRA.....	46
HRV.....	51
HUN.....	56
HYE.....	61
ITA.....	66
KAT.....	71
KAZ.....	76
LAV.....	81
LIT.....	86
NLD.....	91
POL.....	96
POR.....	101
RON.....	106
RUS.....	111
SLK.....	116
SLV.....	121
SPA.....	126
SRP.....	131
UKR.....	136
UZB.....	141

A

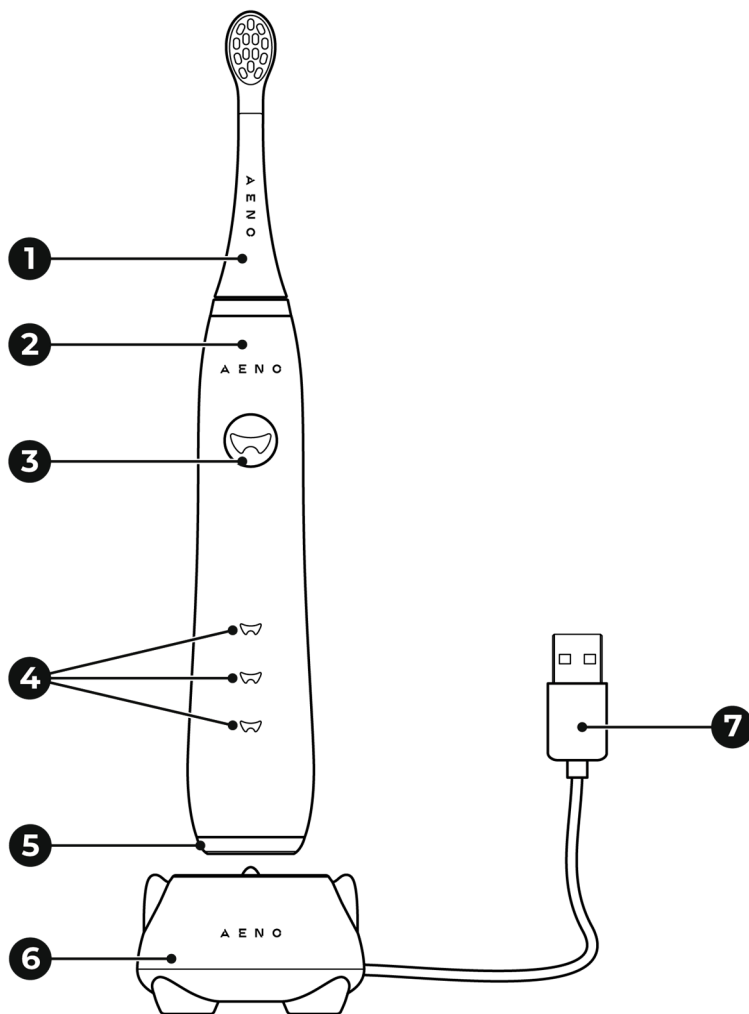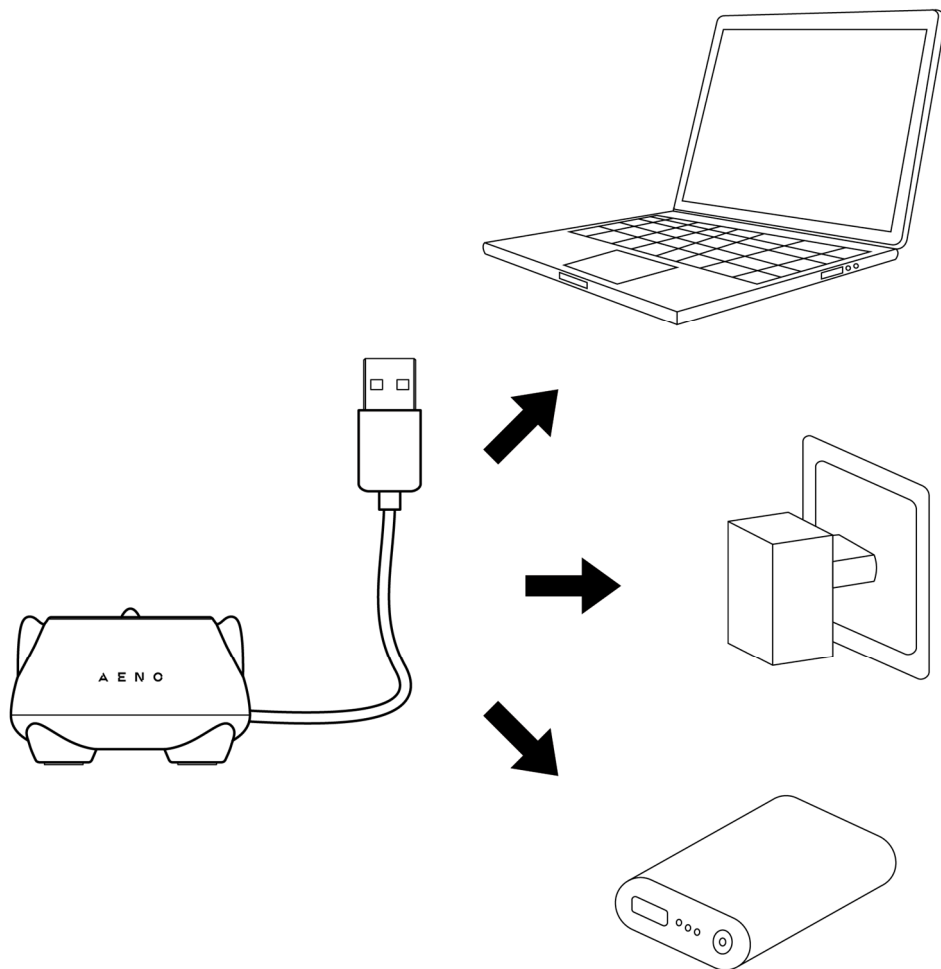
C

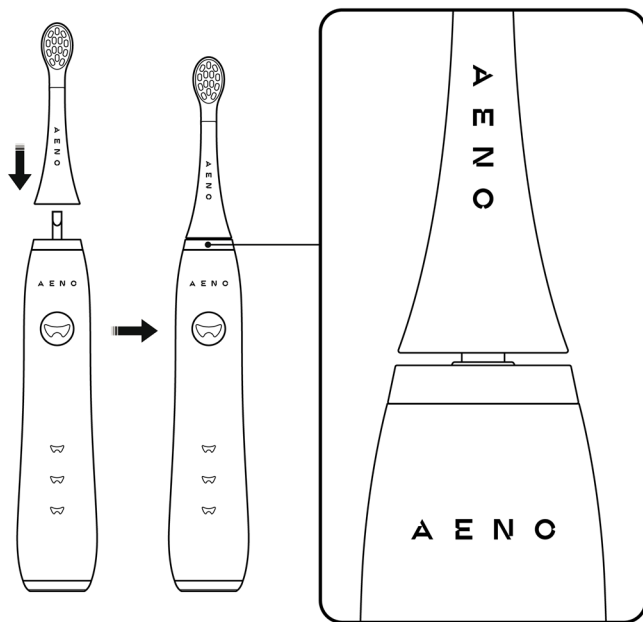D

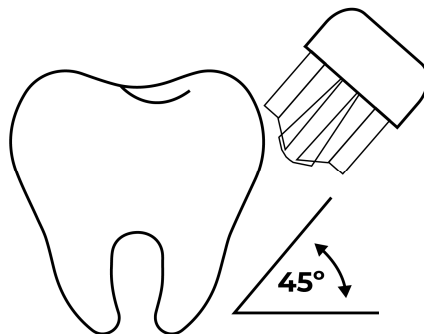ENG **AENO™** Toothbrush is intended for daily oral hygiene.

Model: ADB0009.

Articles: ADB0009V (Violet), ADB0009M (Mint).

Technical Specifications

Toothbrush. Battery: non-removable, Li-ion, 3.7 V, 600 mAh. Full charge time: up to 3 h. Battery life: up to 3 h. Number of operating modes: 3. Vibration intensity: up to 19,200 movements per minute. Built-in pressure sensor. Ingress protection rating: IP67. Noise level (at a distance of 0.3 m): ≤ 65 dB. Control: manual. Toothbrush dimensions without the brush head (diameter × height): 31×137 mm. Weight (without the brush head): 93 g. Casing material: ABS plastic, TPE. Bristle stiffness: soft. Service life of brush heads: 3 months. Operating conditions: temp. +15...+45 °C, RH ≤ 95 %. Storage conditions: temp. +15...+45 °C, RH ≤ 85 %.

Charging Station. Power supply: input 5.0 V / 1.0 A (DC). Ingress protection rating: IP67. Material: ABS plastic. Dimensions (L×W×H): 68×56×37 mm. Weight: 74 g. Connection: USB (A) cable, 1.2 m.

Package Contents

Sonic Toothbrush, charging station with USB (A) cable, replacement brush heads (3 pcs.), Quick Start Guide, game book.

Device Elements and Accessories (see Fig. A)

1 – replacement brush head, 2 – handle, 3 – control button, 4 – indicators, 5 – light ring, 6 – charging station, 7 – USB (A) cable.

Limitations and Warnings

Read this document thoroughly before using the device. Save these instructions for future reference.

If you experience any issues or difficulties while using your AENO™ device, please email support at support@aeno.com or chat with us online at aeno.com/service-and-warranty. The specialists will help you solve the problem, and you will not have to spend time and effort visiting the store.

ATTENTION! If the user has a pacemaker, oral health conditions, or has recently undergone oral surgery, consult a dentist before first use. Consult a pediatric dentist before a child uses the device for the first time.

This toothbrush is a personal hygiene device and is not intended for use on multiple patients in a dental practice or institution. Each family member should use an individual brush head. Replace brush heads at least once every 3 months. Use only original brush heads that are compatible with this model. Use the device only for brushing teeth. The device is not intended for use by children or by persons with physical, sensory, or mental disabilities unless they have sufficient experience or knowledge of how to use the device, or are supervised by a person responsible for their safety. Usage by children under the age of 3 is prohibited. For more information, visit aeno.com/documents

strongly not recommended. Children over the age of 3 should use the device only under the supervision of a person responsible for their safety. Do not allow children to play with the accessories, the toothbrush, or to chew on it or bite it. Do not drop or throw the device. Do not turn on or use a damaged device. Do not expose the device to high temperatures (from heating appliances or direct sunlight). If your gums begin to bleed after using the device and the bleeding continues for several days, contact your dentist. Use only the charging station included in the package to charge the device. Keep the charging station with the cable out of reach of children. Make sure that the voltage and frequency specified in the technical documentation match the parameters of the power source. If you do not use the device for a long time, charge it every 2–3 months to prevent deep discharge of the battery. To extend battery life, do not store the device on a charging station that is connected to a power source. Do not use chemically aggressive detergents, abrasive pastes, products containing acids and solvents, or steel wool sponges to clean the device. Children should not clean the device. Do not wash the device or its accessories in a dishwasher. For a detailed description of the device, its modes and features, see the full User Manual available at aeno.com/documents.

Control Button Functions

Action	Description
Press once (device is turned off)	Turn the device on
Press once (less than 8 seconds have passed since the device was turned on)	Switch operating modes
Press once (more than 8 seconds have passed since the device was turned on)	Turn the device off
Press and hold for 5 seconds (device is turned off)	Deactivate/activate the pressure sensor

Pressure Sensor

The toothbrush has a built-in pressure sensor that warns of excessive pressure on the tooth enamel during brushing. If this happens, the device reduces the vibration intensity and starts vibrating once per second, while the light ring (5) flashes red (see “Indication”).

The pressure sensor is active by default. To deactivate/activate the sensor, the device must be turned off. Press and hold the button (3) for 5 seconds – the light ring will flash in two colors (depending on the pressure sensor status, see “Indication”) for 5 seconds.

Operating Modes

- Soft mode is recommended for beginners and for brushing sensitive teeth;
- Medium mode is recommended for brushing teeth with normal sensitivity;
- Intensive mode is recommended for a more thorough plaque removal from the tooth surface.

Indication

Indicators				
				
The device is working in soft mode	The device is working in medium mode	The device is working in intensive mode	Battery is charging	Battery is fully charged
Light ring				
 R	Excessive pressure on the tooth enamel during brushing			
	The device is turned on, battery level is low			
 R	The battery is fully discharged, the device does not turn on			
 B+G	Pressure sensor is activated	 B+R	Pressure sensor is deactivated	
 C	The device is working in soft mode, pressure sensor is active	 C	The device is working in soft mode, pressure sensor is inactive	
 G	The device is working in medium mode, pressure sensor is active	 G	The device is working in medium mode, pressure sensor is inactive	
 B	The device is working in intensive mode, pressure sensor is active	 B	The device is working in intensive mode, pressure sensor is inactive	

Symbols Description

	Indicator / light ring is on	Light ring color	
	Indicator / light ring is flashing	C	Cyan
	Indicator is off	G	Green
	Light ring flashes for 2 seconds	B	Blue
	Light ring flashes for 5 seconds	R	Red

Preparing for Use

- Carefully open the package and remove the device and the accessories. Remove packing materials.
- Clean the device and accessories (see “Cleaning and Maintenance”).

3. Place the charging station (6) on a dry, flat surface and connect it to a power source (see Fig. **B**).

Note. Power adapter is not included. Use only a safe extra-low voltage (SELV) power adapter that has output parameters of 5.0 V / 1.0 A (DC). The normal recommended temperature for battery charging is +15...+45 °C.

4. Place the handle (2) on the charging station. During charging, the indicators (4) will flash alternately.

When the indicators start glowing steadily, the battery is fully charged.

5. Disconnect the charging station from the power source.

6. Remove the handle from the charging station and attach the brush head (1) as shown in Figure **C**.

Use

1. Wet the bristles of the brush head and apply the desired amount of a dental care product to the bristles.

2. Point the bristles of the brush head toward the tooth surface at an angle of about 45° (see Fig. **D**).

3. Press the button once to turn on the device.

4. Continue to press the button to switch operating modes.

Note. The operating mode can be changed within 8 seconds after the device has been turned on.

5. Brush the teeth without applying excessive pressure to the tooth enamel.

6. Every 30 seconds, the device will notify you to change the brushing zone by briefly stopping the vibration. Change the brushing zone after each notification. After 2 minutes of operation, the device will turn off automatically.

Note. The last selected operating mode is stored in the device memory and will be activated the next time the device is turned on.

Cleaning and Maintenance

ATTENTION! Children should not clean the device.

ATTENTION! Turn off the device and disconnect the charging station from the power source before cleaning the device and accessories.

ATTENTION! Do not wash the device or its accessories in a dishwasher.

Cleaning the Brush Head. Rinse the brush head, including the bristles, under running water.

Cleaning the Handle. Rinse the handle under running water or wipe with a soft cloth dampened with water. Then wipe dry.

Cleaning the charging station. Wipe the station with a soft, dry cloth.

Troubleshooting

The toothbrush does not turn on. Possible cause: the battery is discharged. Solution: charge the device (charging time – up to 3 hours).

The toothbrush is not charging. Possible causes: the charging station is not connected to a power source; the USB (A) cable is damaged. Solution: connect the station to a power source (see Fig. **B**); check the integrity of the USB (A) cable.

Weak vibration of the brush head. Possible cause: the place where the handle and the brush head connect is dirty. Solution: clean the handle and the brush head.

ATTENTION! If none of the suggested solutions helped you eliminate the problem, contact the supplier or an authorized service center. Do not disassemble the device or attempt to repair it yourself.

ASBISC reserves the right to modify the device and to revise and change this document without prior notice to users.

The warranty period and service life are 2 years from the date of retail sale of the device.

Manufacturer: ASBISC Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Tel.: +48732080077, info@asbis.com, asbis.com. Made in China. All trademarks used herein are the property of their respective owners. The date of manufacture is marked on the package.

Up-to-date information about and a detailed description of the device, as well as connection instructions, certificates, and details about the companies that deal with quality and warranty claims are available at aeno.com/documents. This device complies with applicable safety requirements.

RoHS Directive. The device complies with the requirements of the Directive 2011/65/EU on the Restriction of Hazardous Substances (RoHS), including the requirements of the RoHS Directive 2015/863/EU.

Information on disposal. These symbols indicate that when disposing of the device, its batteries and accumulators, as well as its electrical and electronic accessories, you must follow the rules for handling waste electrical and electronic equipment (WEEE) and the rules for handling battery and accumulator waste. According to the regulations, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories with unsorted household waste, as this will harm the environment. This device has a non-removable battery; at the end of the battery's service life, the device must be disposed of. To dispose of this device, return it to the point of purchase or take it to a local recycling center. For detailed information, contact your local household waste disposal service.

For devices with non-removable batteries: under normal or abnormal operating conditions, as well as in the event of a single malfunction, the charging voltage and current must not exceed their maximum permissible values. If the battery temperature exceeds the maximum allowed temperature, the device will stop charging under abnormal conditions. If the battery temperature drops below the minimum allowed temperature during charging, the current will be limited under abnormal conditions.

In case of battery leakage immediately stop using the device. Put on latex/rubber gloves and a protective mask, and open the windows in the room. Place the device with the leaking battery in one or two transparent plastic bags. Dispose of the device in accordance with local recycling regulations.

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