A E N O



Quick Start Guide

Sonic Toothbrush DB7 / DB8


(Models: ADB0007,
ADB0008)

v 2.4.0 | 15.09.2026


ENG	5	KAZ.....	46
ARA	8	LAV	49
BOS	10	LIT	52
BUL.....	13	NLD.....	55
CES	16	POL.....	58
DEU	19	POR.....	61
ELL.....	22	RON	64
EST.....	25	RUS	67
FRA.....	28	SLK.....	70
HRV	31	SLV	73
HUN.....	34	SPA.....	76
HYE.....	37	SRP	79
ITA	40	UKR.....	82
KAT	43	UZB.....	85


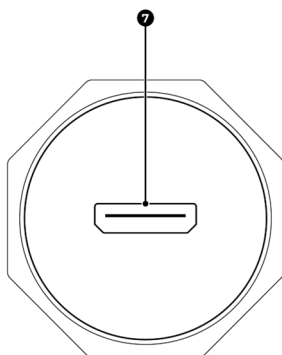
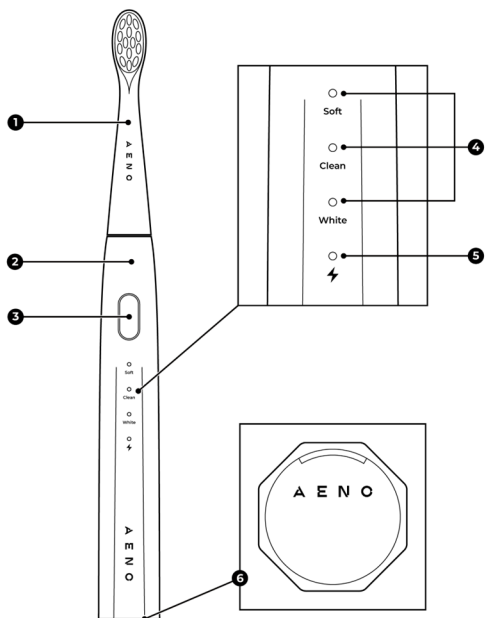
B

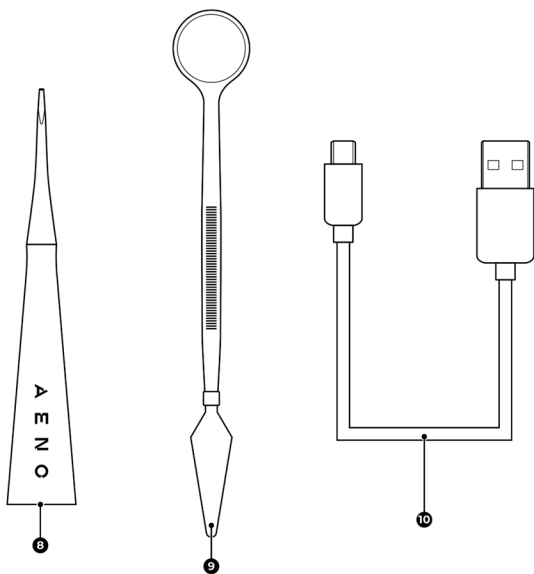C

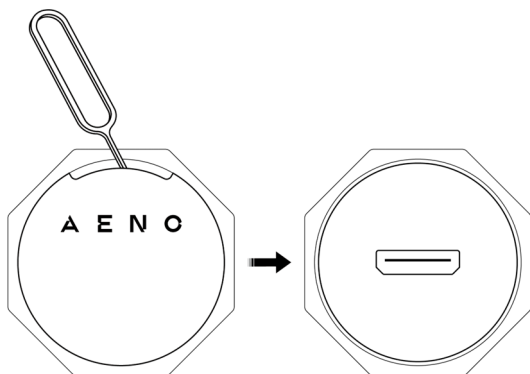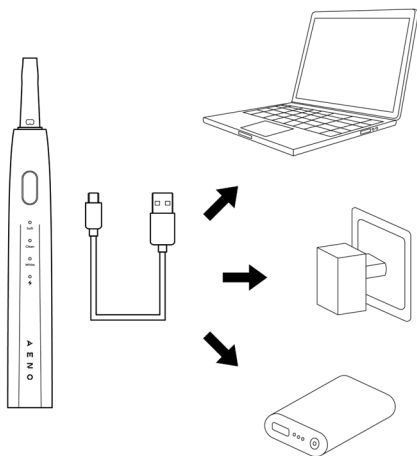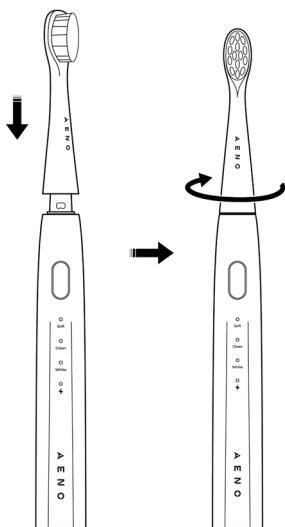D**E**

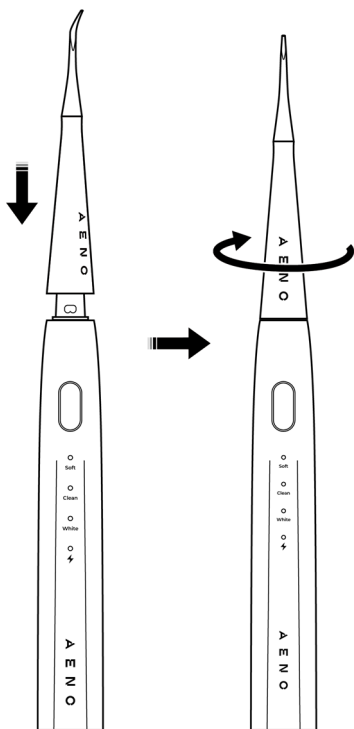ENG **AENO™** Toothbrush is intended for daily oral hygiene.

Models: ADB0007, ADB0008.

Technical Specifications

Power supply: input 5.0 V / 0.3 A (DC). Battery: non-removable, Li-ion, 3.7 V, 800 mAh. Full charge time: up to 3 h. Operating time on a single charge: up to 10 h. Number of operating modes: 3. Vibration intensity: up to 30,000 vibrations per minute. Ingress protection rating: IPX7. Noise level (at a distance of 15 cm): ≥ 55 dB. Control: manual. Color: White. Toothbrush dimensions without the brush head (diameter × height): 21×163 mm. Weight (without the brush head): 49 g. Casing material: ABS plastic, TPE. USB (A) / micro-USB cable length: 30 cm. Bristle stiffness: medium. Service life of brush heads: 3 months. Operating conditions: temp. 0...+40 °C, RH 0–95 %. Storage conditions: temp. 0...+40 °C, RH 0–75 %.

Package Contents (see Fig. A, B)

ADB0007: Sonic Toothbrush, replacement brush head with bristles (1 pc.), USB (A) / micro-USB cable, sticker sheets (2 pcs.), Quick Start Guide.

ADB0008: Sonic Toothbrush, replacement brush heads with bristles (3 pcs.), plaque remover, dental mirror, USB (A) / micro-USB cable, Quick Start Guide.

Device Elements and Accessories (see Fig. A, B)

1 – replacement brush head with bristles, 2 – handle, 3 – control button, 4 – operating mode indicators, 5 – charge indicator, 6 – protective cover, 7 – micro-USB port, 8 – plaque remover, 9 – dental mirror, 10 – USB (A) / micro-USB cable.

Limitations and Warnings

Read this document thoroughly before using the device. Save these instructions for future reference.

If you experience any issues or difficulties while using your AENO™ device, please email support at support@aeno.com or chat with us online at aeno.com/service-and-warranty. The specialists will help you solve the problem, and you will not have to spend time and effort visiting the store.

ATTENTION! If you have a pacemaker, oral health conditions, or have recently undergone oral surgery, consult your dentist before first use.

The device is not intended for use by children or by persons with physical, sensory, or mental disabilities unless they have sufficient experience or knowledge of how to use the device, or are supervised by a person responsible for their safety. The plaque remover is intended for use by adults only. Do not let children play with the device and its accessories. Do not drop or throw the device. Do not turn on or use a damaged device. Do not expose the device to high temperatures (from heating appliances or direct sunlight). If your gums begin to bleed after using the device and the bleeding continues for several days, contact your dentist. Make sure that the power supply parameters specified in the technical specifications match the parameters of the power source. Each family member should use an individual brush head. Replace the brush head with bristles at least once every 3 months. Use only original brush heads designed for these models. If you do not use the device for a long time, charge it every 2–3 months to prevent deep discharge of the battery. Do not use chemically aggressive detergents, abrasive pastes, products containing acids and solvents, or steel wool sponges to clean the device. Do not wash the device and its accessories in a dishwasher. For a detailed description of the device, its modes and features, see the full User Manual available at aeno.com/documents.

Control Button Functions

Action	Description
Press once (device is turned off)	Turn the device on
Press once (less than 5 seconds have passed since the device was turned on)	Switch operating modes
Press once (more than 5 seconds have passed since the device was turned on)	Turn the device off

Statues of the Battery Charge Indicator

Color	Status	Meaning
Green	On	The device is turned on, battery level is medium/high
		The device is connected to a power source; the battery is fully charged
Yellow	On	Battery is charging
Red	On	The device is turned on, battery level is low
	Blinking	The battery is fully discharged, the device does not turn on
–	Off	The device is not connected to a power source or is turned off

Operating Modes Indicators

You can tell the current operating mode of the toothbrush by the steadily lit white indicator located above the name of the mode.

Operating Modes

- “Soft” is recommended for beginners and for brushing sensitive teeth;
- “Clean” is recommended for brushing teeth with normal sensitivity;
- “White” is recommended for teeth whitening. Suitable for people who smoke or frequently consume beverages that affect the color of tooth enamel (tea, coffee, etc.).

Preparing for Use

1. Carefully open the package and remove the device and the accessories. Remove packing materials.
2. Clean the device and accessories (see “Cleaning and Maintenance”).
3. Remove the protective cover (6) from the toothbrush handle (see Fig. C).

ATTENTION! The protective cover ensures that the micro-USB port (7) is waterproof. If the protective cover is removed, do not allow liquids to enter the device port.

Note. The power adapter is not included. To charge the device, use a power adapter with output parameters of 5.0 V / 1.0 A (DC).

4. Connect the device to the power source using the USB (A) / micro-USB cable (10) (see Fig. D). The charge indicator (5) will light up yellow. When the indicator turns green, the battery is fully charged.
5. Disconnect the device from the power source, then disconnect the cable from the micro-USB port. Insert the protective cover into the device handle.

Use

1. Select the required brush head:
 - **Brush head with bristles.** Install the brush head (1) as shown in Fig. E. Wet the bristles of the brush head and apply the desired amount of a dental care product to the bristles.
 - **Plaque remover.** Install the plaque remover (8) as shown in Fig. F.
2. Direct the bristles / metal tip of the brush head toward the tooth surface at an angle of approximately 45°.
3. Press the button (3) once to turn on the device.
4. Continue to press the button to switch operating modes.

Note. The operating mode can be changed within 5 seconds after the device has been turned on.

5. Brush the teeth without applying excessive pressure to the tooth enamel. When using the plaque remover, clean problem areas with gentle vertical movements.

ATTENTION! Do not apply excessive force when using the plaque remover. Careless use of the remover can cause oral trauma as well as damage tooth enamel.

6. Every 30 seconds, the device will notify you to change the brushing zone by briefly stopping the vibration. Change the brushing area after each notification. After 2 minutes of operation, the device will turn off automatically.

Note. The last selected operating mode is stored in the device memory and will be activated the next time the device is turned on.

Note. Using a plaque remover cannot substitute professional oral hygiene at the dentist.

Cleaning and Maintenance

ATTENTION! Before cleaning the device and accessories, turn off the device and disconnect it from the power source. Make sure that the protective cover is securely inserted into the handle.

ATTENTION! Do not wash the device and its accessories in a dishwasher.

Rinse the device and accessories under running water. Then wipe dry.

Troubleshooting

The toothbrush does not turn on. Possible cause: battery is discharged. Solution: charge the device (charging time – up to 3 hours).

The toothbrush is not charging. Possible causes: the power source is faulty; the USB (A) / micro-USB cable is damaged. Solution: connect the device to another power source (see Fig. D); check the USB (A) / micro-USB cable for damage.

Weak vibration of the brush head. Possible cause: the place where the handle and the brush head connect is dirty. Solution: clean the handle and the brush head.

ATTENTION! If none of the suggested solutions helped you eliminate the problem, contact the supplier or an authorized service center. Do not disassemble the device or attempt to repair it yourself.

Information contained herein may be changed without prior notice to users. For the latest information about the device, certificates, warranty, and quality, see the complete User Manual at aeno.com/documents. This device complies with applicable safety requirements.

Warranty. The warranty period and service life of the device are 2 years from the date of retail sale of the device, unless otherwise required by local law. Additional warranty information is available at aeno.com/service-and-warranty.

Manufacturer: ASBISC Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Tel.: +48732080077, info@asbis.com, asbis.com. Made in China. The date of manufacture is indicated on the package. All trademarks and trade names used herein are the property of their respective owners.

RoHS Directive. The device complies with the requirements of the Directive 2011/65/EU on the Restriction of Hazardous Substances (RoHS), including the requirements of the RoHS Directive 2015/863/EU.

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ENG Disposal information. These symbols indicate that when disposing of the device, its batteries and accumulators, as well as its electrical and electronic accessories, you must follow the rules for handling waste electrical and electronic equipment (WEEE) and the rules for handling battery and accumulator waste. According to the regulations, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories with unsorted household waste, as this will harm the environment. This device has a non-removable battery; at the end of the battery's service life, the device must be disposed of. To dispose of this device, return it to the point of purchase or take it to a local recycling center. For detailed information, contact your local household waste disposal service.

For devices with non-removable batteries: under normal or abnormal operating conditions, as well as in the event of a single malfunction, the charging voltage and current must not exceed their maximum permissible values. If the battery temperature exceeds the maximum allowed temperature, the device will stop charging under abnormal conditions. If the battery temperature drops below the minimum allowed temperature during charging, the current will be limited under abnormal conditions.

In case of battery deformation or leakage, immediately stop using the device. Put on latex/rubber gloves and a protective mask and open the windows in the room. Place the device with the leaking battery in one or two transparent plastic bags. Dispose of the device in accordance with local recycling regulations.

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