
A E N O

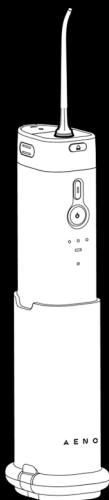Quick Start Guide



Dental Irrigator DI2

(Model: ADI0002)


v 3.2.3 | 21.08.2026

ENG	4
ARA	7
BOS	9
BUL.....	12
CES	16
DEU	19
ELL.....	23
EST.....	26
FRA.....	29
HRV	33
HUN.....	36
HYE.....	40
ITA	43
KAT.....	47


KAZ	51
LAV.....	55
LIT.....	58
NLD	61
POL	65
POR.....	69
RON.....	73
RUS	77
SLK	81
SLV	84
SPA.....	87
SRP.....	91
UKR.....	94
UZB	98

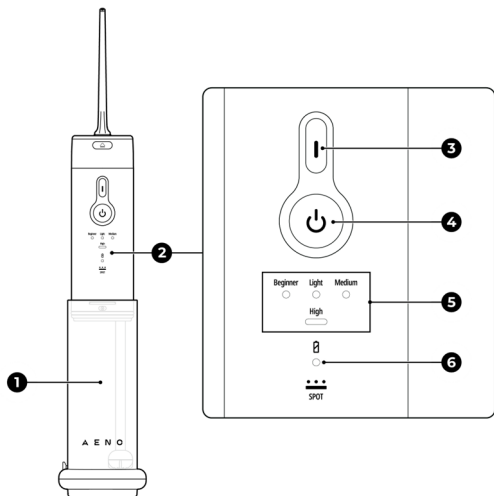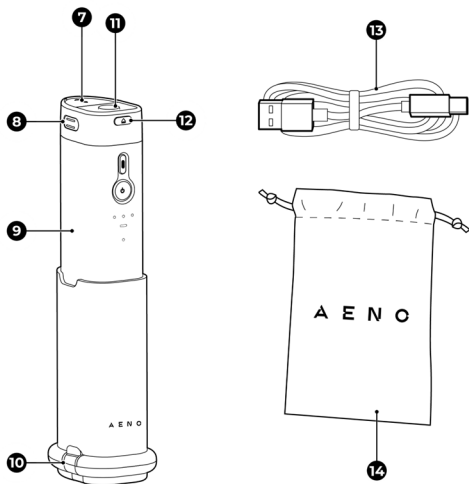A**B**

C

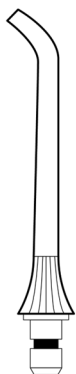15

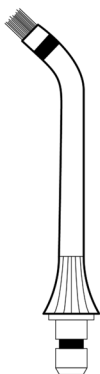16




17

D

E

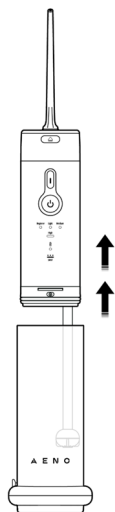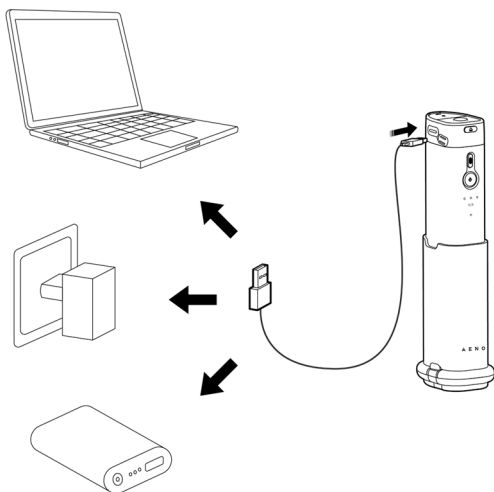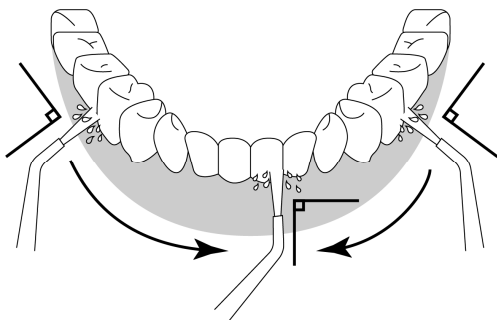F**G**

ENG The **ADI0002** Dental Irrigator is used to remove dental plaque and food debris from hard-to-reach areas of the mouth.

Technical Specifications

Rated voltage/power: 3.7 V / 8.0 W (DC). Charging voltage/current (input): 5.0 V / 2.0 A (max.) (DC). Battery: lithium-ion, non-removable, 3.7 V, 2,500 mAh, 9.25 Wh. Water pressure: 413.7–965.3 kPa. Water pulsation: 1,900 pulsations/min. Water reservoir: 200 ml. Full battery charge time: up to 3.5 h. Operating time on a single charge: at least 200 min. Operating time on a full reservoir: at least 40 s. Ingress protection rating: IPX7. Noise level (at a distance of 0.3 m): ≤ 72 dB. Control: manual. Color: White. USB (A) / USB (C) cable length: 0.6 m. Material: ABS plastic. Dimensions, collapsed (L×W×H): 73×53×136 mm. Weight (without nozzles): 253 g. Operating conditions: temp. +5...+40 °C, RH 30–75 %. Storage conditions: temp. –10...+60 °C, RH ≤ 75 %.

Package Contents (see Fig. B, C)

ADI0002 Dental Irrigator, removable nozzles (3 pcs.), USB (A) / USB (C) cable, travel bag, Quick Start Guide.

Device Elements and Accessories (see Fig. A, B, C)

1 – water reservoir, 2 – control panel, 3 – operating mode selection button, 4 – power button, 5 – operating mode indicator, 6 – battery charge indicator, 7 – water reservoir fill door, 8 – USB (C) port with a cover, 9 – main body, 10 – water drain hole with a plug, 11 – opening for the nozzle, 12 – nozzle release button, 13 – USB (A) / USB (C) cable, 14 – travel bag, 15 – standard nozzle, 16 – orthodontic nozzle, 17 – periodontal nozzle.

Limitations and Warnings

The device is not intended for use by children or by persons with physical, sensory, or mental disabilities unless they have sufficient experience or knowledge of how to use the device, or are supervised by a person responsible for their safety. Do not let children play with the device and accessories. If you have oral cavity diseases, you should consult your dentist before using the device. Do not use a damaged device or accessories. If the device has been transported or stored at low temperatures, leave it at room temperature for 2 hours before use. Do not drop or throw the device. Clean the device only after turning it off and disconnecting it from the power supply. Do not use chemical or corrosive cleaning agents, abrasive pastes, agents that contain acids or solvents, or metal sponges to clean the device. Do not wash the device or its accessories in a dishwasher. Only use original accessories. Replace nozzles at least once every 3 months. Do not point the outlet of the nozzle into your eyes while the device is in use. If you are not using the device for a long time, charge it every 2–3 months to prevent deep discharge of the battery. For a detailed description of the device, its modes and features, see the full User Manual available at aeno.com/documents.

Intended Purpose of Nozzles

- Standard nozzle (15) is designed to remove plaque and food debris from between the teeth;
- Orthodontic nozzle (16) is designed to clean dental braces;
- Periodontal nozzle (17) is designed to clean and rinse periodontal pockets.

If you experience any issues or difficulties while using your AENO™ device, please email support at support@aeno.com or chat with us at aeno.com/service-and-warranty. Our specialists will help you solve the problem, and you will not have to spend time and effort visiting the store.

Button Functions

Name	Action	Function
Operating mode selection button (3)	Press once	Select operating mode
Power button (4)	Press once	Turn the device on/off
Nozzle release button (12)	Press and hold	Release the nozzle to remove it

Battery Charge Indication

Color	Status	Description
Green	On	The device is connected to a power supply and fully charged
		The device is turned on. The charge level is > 10 %
Red	Blinking	The device is connected to a power supply. The device is charging
		The device is turned on. The charge level is < 10 %
		The device's battery is fully discharged, the device does not turn on
–	Off	The device is turned off or not connected to a power supply

Operating Mode Indication

When the device is running in "Beginner", "Light", "Medium" or "High" modes, the indicator located on the control panel (2) under the name of the selected mode is solid white.

When the device is running in "Spot" mode, all 4 operating mode indicators on the control panel are solid white.

Operating Modes of the Device

- "Beginner": low water pressure. Suitable for cleaning sensitive teeth and gums, as well as for first-time users.
- "Light": moderate water pressure. Suitable for gentle cleaning. Recommended as the next mode after "Beginner" for a smooth transition to modes with a stronger water pressure.
- "Medium": standard water pressure. Suitable for daily cleaning.
- "High": high water pressure. Suitable for intensive cleaning. Recommended for experienced users.
- "Spot": thorough cleaning of the oral cavity, removal of food debris from interdental spaces and prevention of tartar formation.

Note. In "Spot" mode, the device alternately ejects a jet of water for 2 seconds and then stops the water ejection for 2 seconds.

Preparing the Device for Use

1. Carefully open the package and take out the device and its accessories. Remove the device from the travel bag (14).
2. Pull the main body (9) all the way up (see Fig. **D**).
3. Detach the reservoir (1) as shown in Fig. **E**. Rinse it and the nozzles under running water, then wipe them with a soft, dry cloth.

Note. The reservoir locks in place tightly. You need to apply force to detach it.

4. Reattach the reservoir and put the device on a flat, stable surface.
5. Remove the cover from the USB (C) port (8) and connect the device to a power supply using the USB (A) / USB (C) cable (13) (see Fig. **F**).

Note. Power adapter is not included. Use only a safe extra-low voltage (SELV) power adapter that has output parameters of 5.0 V / 2.0 A (max.) (DC).

6. Fully charge the device (takes at least 3.5 hours). Once the battery is fully charged, the battery charge indicator (6) will turn green.
7. Disconnect the device from the power supply and insert the cover into the USB (C) port of the device.

Use

1. Open the fill door (7) and fill the reservoir with water.

ATTENTION! Only use purified water without adding mouthwashes and other liquids or substances.

2. Close the fill door.
3. Choose the desired nozzle (see "Intended Purpose of Nozzles") and insert it into the opening (11) until it clicks into place.

Note. The nozzles can rotate 360°.

4. Set the device to the desired mode using the operating mode selection button (see "Operating Modes of the Device").
5. Point the outlet of the nozzle into the oral cavity and turn the device on by pressing the power button.
6. Carry out oral hygiene by moving the outlet of the nozzle along the edge of the gum as shown in Figure G.

Note. If the reservoir runs out of water during use, turn off the device by pressing the power button once and refill the reservoir as described in Steps 1 and 2 above. After that, continue performing your oral hygiene.

7. Turn off the device by pressing the power button once.
8. Remove the plug (10) and drain the remaining water.
9. Clean the device (see "Daily care").

Cleaning and Maintenance

Daily care. Direct the outlet of the nozzle from yourself and turn the device on by pressing the power button once. Wait for the device to stop ejecting water completely. After that, turn the device off by pressing the power button again. Press the nozzle release button and pull the nozzle up to remove it. Detach the reservoir as shown in Figure E. Rinse the nozzle and reservoir under running water, then dry them thoroughly. Reattach the reservoir and nozzle.

Total care. The reservoir is to be cleaned thoroughly once a week. Detach the reservoir as shown in Figure E. Wash it outside and inside under running water with a soft brush or sponge, without using cleaning agents, then dry it thoroughly. Reattach the reservoir.

Troubleshooting

The device does not turn on. Possible cause: the battery is discharged. Solution: charge the device's battery (charging time is 3.5 hours minimum).

The device does not eject water. Possible causes: the reservoir is out of water; the outlet of the nozzle is clogged. Solution: fill the reservoir with water; clean the outlet of the nozzle.

Insufficient water pressure. Possible cause: the wrong operating mode has been selected. Solution: set the required operation mode using the mode selection button.

ATTENTION! If none of the suggested solutions helped you eliminate the problem, contact the supplier or an authorized service center. Do not disassemble the device or attempt to repair it yourself.

Information contained herein may be changed without prior notice to users. For the latest information about the device, certificates, warranty, and quality, see the complete User Manual at aeno.com/documents. This device complies with applicable safety requirements.

Warranty. The warranty period and service life of the device are 2 years from the date of retail sale of the device, unless otherwise required by local law. Additional warranty information is available at aeno.com/service-and-warranty.

Manufacturer: ASBISC Enterprises PLC, lapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Tel.: +48732080077, info@asbis.com, asbis.com. Made in China. The date of manufacture is indicated on the package. All trademarks and trade names used herein are the property of their respective owners.

RoHS Directive. The device complies with the requirements of the Directive 2011/65/EU on the Restriction of Hazardous Substances (RoHS), including the requirements of the RoHS Directive 2015/863/EU.

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aeno.com/documents

ENG Disposal information. These symbols indicate that when disposing of the device, its batteries and accumulators, as well as its electrical and electronic accessories, you must follow the rules for handling waste electrical and electronic equipment (WEEE) and the rules for handling battery and accumulator waste. According to the regulations, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories with unsorted household waste, as this will harm the environment. This

device has a non-removable battery; at the end of the battery's service life, the device must be disposed of. To dispose of this device, return it to the point of purchase or take it to a local recycling center. For detailed information, contact your local household waste disposal service.

For devices with non-removable batteries: under normal or abnormal operating conditions, as well as in the event of a single malfunction, the charging voltage and current must not exceed their maximum permissible values. If the battery temperature exceeds the maximum allowed temperature, the device will stop charging under abnormal conditions. If the battery temperature drops below the minimum allowed temperature during charging, the current will be limited under abnormal conditions.

In case of battery deformation or leakage, immediately stop using the device. Put on latex/rubber gloves and a protective mask and open the windows in the room. Place the device with the leaking battery in one or two transparent plastic bags. Dispose of the device in accordance with local recycling regulations.

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