A E N O

Quick Start Guide



Dental Irrigator DI3


(Model: ADI0003)

v 1.1.0 | 07.05.2026


ENG.....	10
ARA.....	16
BOS.....	20
BUL.....	26
CES.....	33
DEU.....	39
ELL.....	46
EST.....	52
FRA.....	58
HRV.....	65
HUN.....	71
HYE.....	78
ITA.....	84
KAT.....	90
KAZ.....	97
LAV.....	104
LIT.....	110
NLD.....	116
POL.....	123
POR.....	130
RON.....	136
RUS.....	143
SLK.....	150
SLV.....	156
SPA.....	162
SRP.....	169
UKR.....	175
UZB.....	182

A

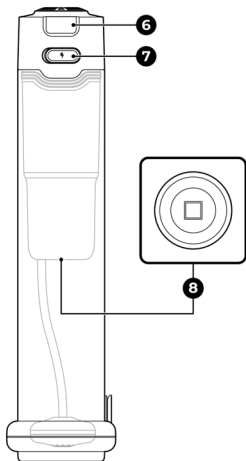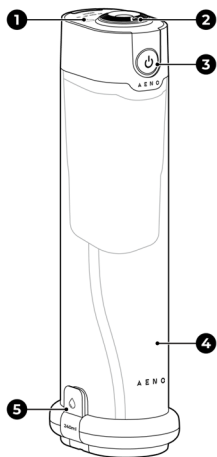
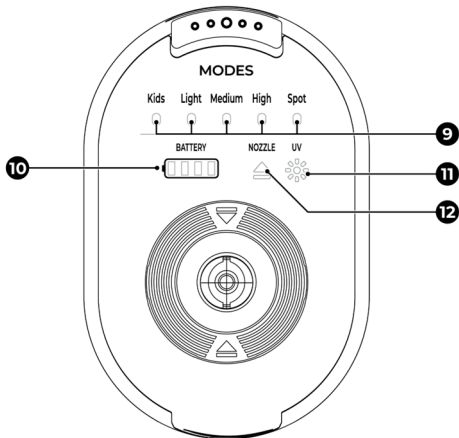
C

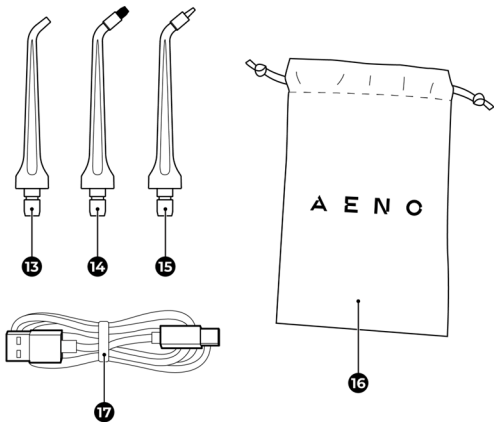
D

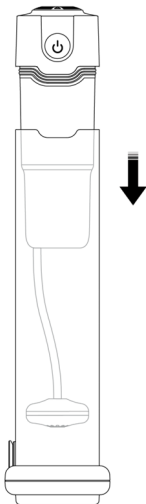E

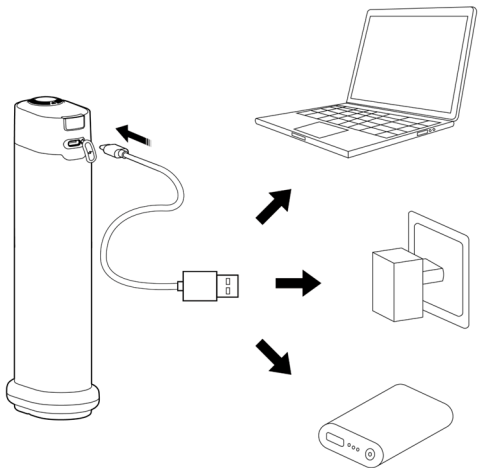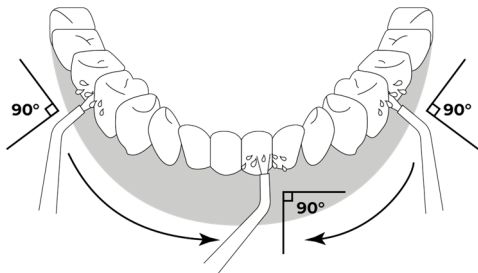ENG The **ADI0003** Dental Irrigator is designed to remove dental plaque and food debris from hard-to-reach areas of the mouth.

Technical Specifications

Power supply (input): 5.0 V / 2.0 A (DC), 10.0 W (max.). Battery: non-removable, lithium-ion, 2,500 mAh, 3.7 V. Water pressure: up to 965 kPa. Pulsation: up to 1,700 pulsations/min. Water tank: removable, 240 ml. UV LED. Full battery charge time: up to 3.5 h. Operating time on a single battery charge: up to 200 min. Operating time on a full tank: at least 55 s. Ingress protection rating: IPX7. Noise level (at a distance of 0.3 m): ≤ 72 dB. Control: manual. Color: White, Black. USB cable length: 0.6 m. Material: ABS plastic. Dimensions without nozzle (L×W×H): 77×54×225 mm. Net weight: 370 g. Operating conditions: temp. +5...+40 °C, RH ≤ 75 %. Storage conditions: temp. 0...+45 °C, RH ≤ 75 %.

Package Contents (see Fig. A, C)

ADI0003 Dental Irrigator, replacement nozzles (3 pcs.), USB (A) / USB (C) cable, travel bag, Quick Start Guide.

Device Elements and Accessories (see Fig. A, B, C)

1 – display, 2 – opening for the nozzle, 3 – power button, 4 – water tank, 5 – water tank plug, 6 – operating mode selection button, 7 – USB (C) port with a cover, 8 – UV LED, 9 – operating mode indicators, 10 – battery charge indicator, 11 – “UV” indicator, 12 – “NOZZLE” indicator, 13 – standard nozzle, 14 – orthodontic nozzle, 15 – periodontal nozzle, 16 – travel bag, 17 – USB (A) / USB (C) cable.

Limitations and Warnings

Read this document thoroughly before using the device. Save these instructions for future reference.

If you experience any issues or difficulties while using your AENO™ device, please email support at **support@aeno.com** or chat with us at **aeno.com/service-and-warranty**. The specialists will help you solve the problem, and you will not have to waste time and effort visiting the store.

ATTENTION! If you have a pacemaker, oral diseases, or if you have just undergone oral surgery, consult your doctor before using the device.

ATTENTION! Do not direct the UV LED rays into the eyes.

ATTENTION! It is forbidden to activate the “UV” function if the water tank is not attached. Make sure that the tank is attached properly before activating the function. Be sure to turn off the “UV” function before detaching the tank.

The device is not intended for use by children or by persons with physical, sensory, or mental disabilities unless they have sufficient experience or knowledge of how to use the device, or are supervised by a person responsible for their safety. Do not allow children to play with the device and its accessories. Do not use a damaged device or accessories. Do not drop or throw the device. Make sure that the voltage and frequency specified in the technical documentation correspond to the power supply parameters. Do not point the outlet of the nozzle towards your eyes while the device is running. Each family member should use an individual nozzle. Replace the nozzles at least once every 3 months. Use only original replacement nozzles. If you are not using the device for a long time, charge it every 2–3 months to prevent deep discharge of the battery. For a detailed description of the device, its modes and features, see the complete User Manual available at aeno.com/documents.

Intended Purpose of Nozzles

- Standard nozzle (13) is designed to remove dental plaque and food debris.
- Orthodontic nozzle (14) is designed to clean dental braces.
- Periodontal nozzle (15) is designed to clean and rinse periodontal pockets.

Button functions

Action	Description
Power button (3)	
Press once	Turn the device on/off
Press and hold	Turn the “UV” function on/off
Operating mode selection button (6)	
Press once	Select operating mode

Display Indicators (see Fig. B)

	On (white)	The device is connected to a power supply and fully charged
		The device is on, charge level is $\geq 10\%$
	Flashing white	The device is charging, charge level is $\geq 10\%$
	Flashing red	The device is on, charge level is $< 10\%$
		The device is charging, charge level is $< 10\%$
	Off	The device is turned off or not connected to a power supply
	Flashing white	The nozzle is not inserted into the opening (2), the device does not turn on
	On (blue)	The "UV" function is enabled
	Off	The "UV" function is disabled
	On (white)	The corresponding operating mode is selected
	Flashing white	The water flow is obstructed, the operation of the device is paused
	Off	The device is turned off

* Schematic view of any indicator (9).

Operating Modes of the Device

- "Kids" is recommended for use by children and first-time users.
- "Light" is recommended for cleaning sensitive teeth and gums.
- "Medium" is recommended for daily oral care.
- "High" is recommended for advanced users.
- "Spot" is recommended for cleaning precise spots of the oral cavity.

Note. In "Spot" mode, the device ejects water for 2 seconds and then stops the water ejection for 2 seconds.

Preparing the Device for Use

1. Carefully open the packaging and remove the device and its accessories. Remove any packaging materials.

ATTENTION! Before detaching the tank, make sure that the “UV” function is turned off (the LED (8) is off).

2. Pull the tank (4) down to detach it (see Fig. **D**). Rinse it and the nozzles under running water, then wipe them with a soft, dry cloth.
3. Reattach the tank and place the device on a flat, stable surface.
4. Open the cover of the USB (C) port (7) and connect the device to a power source using the USB (A) / USB (C) cable (17) (see Fig. **E**).

Note. Power adapter is not included. Use only a safe extra-low voltage (SELV) power adapter that has output parameters of 5.0 V / 2.0 A (DC). The normal recommended temperature for battery charging is +5...+40 °C.

5. Fully charge the device (up to 3.5 hours). When the battery is fully charged, the battery charge indicator (10) will turn solid white.
6. Disconnect the device from the power supply and insert the cover into the USB (C) port of the device.

Use

ATTENTION! Before detaching the tank, make sure that the “UV” function is turned off (the LED is off).

1. Pull the tank down to detach it (see Fig. **D**). Fill the tank with water without exceeding the “MAX” mark.

ATTENTION! Only use purified water without adding mouthwashes and other liquids or substances.

2. Reattach the tank. Choose the desired nozzle (see “Intended Purpose of Nozzles”) and insert it into the opening until it clicks into place.
3. If necessary, activate the “UV” function to disinfect the water in the tank (see “UV Disinfection”).

4. Select the desired mode using the operating mode selection button (see “Operating Modes of the Device”).
5. Direct the outlet of the nozzle into the oral cavity and turn the device on by pressing the power button.
6. Carry out oral hygiene by moving the outlet of the nozzle along the edge of the gum as shown in Figure F.

Note. If the tank runs out of water while in use, turn off the device by pressing the power button once and make sure that the “UV” function is turned off (the LED is off). After, detach and refill the tank as described in step 1. Then resume your oral hygiene.

7. Turn off the device by pressing the power button once.
8. Make sure that the “UV” function is turned off (the LED is off). Open the plug (5) and drain the remaining water.
9. Clean the device (see “Daily care”).

UV Disinfection

This function allows you to disinfect the tank and the water in it using UV rays. Press and hold the power button for 2 seconds to turn on the UV LED. The indicator (11) will turn blue. The “UV” function can be activated while the device is on or off. The LED will turn off automatically within 5 minutes of being switched on. To turn it off manually, press and hold the power button for 2 seconds. For more information about this function, refer to the complete User Manual available at aeno.com/documents.

Cleaning and Maintenance

ATTENTION! Before cleaning the device, be sure to turn it off or disconnect it from the power supply and check that the “UV” function is turned off (the LED is off).

ATTENTION! Do not wash the device or its accessories in a dishwasher.

ATTENTION! Do not use chemical and aggressive detergents, abrasive pastes, products containing acids and solvents, or metal sponges to clean the device and accessories.

Daily care. Point the outlet of the nozzle into the sink, then turn the dental irrigator on by pressing the power button once. Wait for the device to stop ejecting water completely. After that, turn the device off by pressing the power button again. Remove the nozzle by pulling it up, and then detach the tank (see Fig. D). Rinse the nozzle and the tank under running water, then dry them thoroughly. Reattach the tank and the nozzle. If necessary, disinfect the tank using UV rays (see “UV Disinfection”).

Total care. Clean the tank thoroughly once a week. Detach the tank (see Fig. D). Wash it inside and out under running water with a soft brush or sponge, without the use of detergents or cleaning agents, and then dry it thoroughly. Reattach the tank. If necessary, disinfect the tank using UV rays (see “UV Disinfection”).

Troubleshooting

The device does not turn on. Possible cause: the battery is discharged. Solution: charge the device (charging time – at least 3.5 hours).

The device does not eject water. Possible causes: the tank is out of water; the outlet of the nozzle is clogged. Solution: fill the tank with water; clean the outlet of the nozzle.

Insufficient water pressure. Possible cause: the wrong operating mode has been selected. Solution: select the desired operating mode using the mode selection button.

ATTENTION! If none of the suggested solutions helped you eliminate the problem, contact the supplier or an authorized service center. Do not disassemble the device or attempt to repair it yourself. Any repairs related to UV LED should be carried out by a qualified technician only using appropriate protective equipment.

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ENG ASBISC reserves the right to modify the device and to revise and change this document without prior notice to users.

The warranty period and service life are 2 years from the date of retail sale of the device.

Manufacturer: ASBISC Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Tel.: +48732080077, info@asbis.com, asbis.com. Made in China. All trademarks used herein are the property of their respective owners. The date of manufacture is marked on the package.

Up-to-date information about and a detailed description of the device, as well as connection instructions, certificates, and details about the companies that deal with quality and warranty claims are available at **aeno.com/documents**. This device complies with applicable safety requirements.

RoHS Directive. The device complies with the requirements of the Directive 2011/65/EU on the Restriction of Hazardous Substances (RoHS), including the requirements of the RoHS Directive 2015/863/EU.

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ENG Information on disposal. These symbols indicate that when disposing of the device, its batteries and accumulators, as well as its electrical and electronic accessories, you must follow the rules for handling waste electrical and electronic equipment (WEEE) and the rules for handling battery and accumulator waste. According to the regulations, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories with unsorted household waste, as this will harm the environment. This device has a non-removable battery; at the end of the battery's

service life, the device must be disposed of. To dispose of this device, return it to the point of purchase or take it to a local recycling center. For detailed information, contact your local household waste disposal service.

For devices with non-removable batteries: under normal or abnormal operating conditions, as well as in the event of a single malfunction, the charging voltage and current must not exceed their maximum permissible values. If the battery temperature exceeds the maximum allowed temperature, the device will stop charging under abnormal conditions. If the battery temperature drops below the minimum allowed temperature during charging, the current will be limited under abnormal conditions.

In case of battery leakage, immediately stop using the device. Put on latex/rubber gloves and a protective mask and open the windows in the room. Place the device with the leaking battery in one or two transparent plastic bags. Dispose of the device in accordance with local recycling regulations.

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